

GA-A22340

CONF 960543--42

**PLASMA MASS DENSITY, SPECIES MIX  
AND FLUCTUATION DIAGNOSTICS RECEIVED  
USING FAST ALFVÉN WAVE**

FEB 0 3 1997

OS

by  
H. IKEZI, J.S. deGRASSIE, R.I. PINSKER,  
and R.T. SNIDER

JUNE 1996

**MASTER**

 **GENERAL ATOMICS**

## DISCLAIMER

This report was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor any agency thereof, nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency thereof. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof.

**DISCLAIMER**

**Portions of this document may be illegible  
in electronic image products. Images are  
produced from the best available original  
document.**

GA-A22340

**PLASMA MASS DENSITY, SPECIES MIX  
AND FLUCTUATION DIAGNOSTICS  
USING FAST ALFVÉN WAVE**

by  
**H. IKEZI, J.S. deGRASSIE, R.I. PINSKER,  
and R.T. SNIDER**

This is a preprint of a paper to be presented at the 11th  
Topical High Temperature Plasma Diagnostics  
Conference, May 12-16, 1996, in Monterey, California,  
and to be published in *Rev. Sci. Instrum.*

**Work supported by  
the U.S. Department of Energy  
under Contract No. DE-AC03-89ER51114**

**GA PROJECT 3466  
JUNE 1996**

DISTRIBUTION OF THIS DOCUMENT IS UNLIMITED



*RF*

## ABSTRACT

We propose to employ a fast Alfvén wave interferometer and reflectometer as a tokamak diagnostic to measure the plasma mass density, D-T species mix profile, and density fluctuations. Utilize the property that the phase velocity of the fast wave propagating across the magnetic field is the Alfvén speed with thermal correction, this fast wave interferometer on the DIII-D tokamak was successfully used to obtain the line integrated density. Since the position of the ion-ion hybrid cut-off in tokamaks is uniquely determined by the species mix ratio and the wave frequency, the reflectometer arrangement finds the species mix profile. The inversion method of reflectometry is discussed. The multiple chord interferometer also measures the mass density fluctuation profile.



## I. MASS DENSITY

The phase velocity of the fast wave propagating across the magnetic field equals the Alfvén speed,  $v_A = B/(4\pi\rho)^{1/2}$ , with the plasma pressure correction

$$\omega/k_{\perp} = v_A (1 + \beta)^{1/2} \quad (1)$$

if the frequency is higher than a few times the ion cyclotron frequencies,  $\omega/\Omega_i \geq 3$ . The measurement of wave phase (wavelength) by interferometry gives the mass density,  $\rho$ . The appropriate frequency for a D-T plasma of density =  $10^{14}/\text{cc}$  with  $B = 4.5$  T is typically 200 MHz which gives the wavelength in the plasma of 3 cm.

The wave sensor can be a single turn  $\dot{B}$  loop of 1–2 cm in diameter. A loop antenna can also be used for the wave launcher. The antenna and the sensor loops can be placed behind the first wall with the signals propagating through the gaps in the wall. The required power level is a few watts. The density profile data may be obtained by employing an array of sensor loops.

The preliminary tests of the fast wave interferometer have been carried out on DIII-D and indicate good agreement with a CO<sub>2</sub> laser interferometer.<sup>1</sup> An example of the results is shown in Fig. 1. The

difference of the results from two interferometers at higher density, where the neutral beam is injected, is partially due to the fact that the cord locations are different. The parameters employed are:  $\omega/2\pi = 60$  MHz,  $B = 1$  T and 2 T, transmitter power = 10 W. The  $\dot{B}$  loops used survive up to 500°C.

The frequency can be swept in a wide range very easily at the fast wave frequency range. The frequency dependence of the phase shift determines the density instantaneously. This operation of interferometer avoids the fringe counting errors which occur in the laser interferometers.

This is a simple diagnostic which does not have the vibration elimination difficulties arising in CO<sub>2</sub> laser system. Since the frequency is in the VHF band and the sensors and antenna are only  $\dot{B}$  loops, this is a low cost diagnostic which can give real time response.



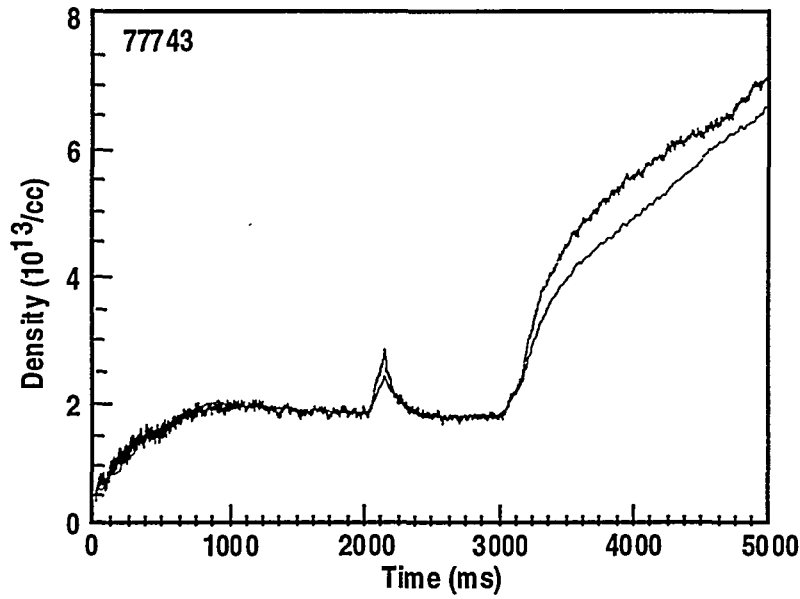


Fig. 1. Densities measured by both the CO<sub>2</sub> interferometer (noisier curve) and by the Alfvén speed (quieter curve).



## II. D-T CONCENTRATION RATIO

The fast wave has ion-ion hybrid resonance(s) if the plasma consists of two or more ion species. The resonance and cut-off frequencies are uniquely determined by the species mix ratio.

The index of refraction perpendicular to B is given by  $N_{\perp}^2 = RL/S$  in Stix's notation.<sup>2</sup> The resonance and cut-off frequencies, determined by  $S = 0$  and  $R = L = 0$ , are

$$\begin{aligned}\omega_{\text{res}} &= \left[ \frac{3-f}{3+(3/2)f} \right]^{1/2} \Omega_D \\ \omega_{\text{cut}} &= \left( 1 - \frac{1}{3}f \right) \Omega_D\end{aligned}\tag{2}$$

in a D-T plasma with the deuterium concentration fraction,  $f$ . Here,  $\Omega_D$  is the ion cyclotron frequency of deuterium. In a non-uniform magnetic field,  $B \propto 1/R$ , in a tokamak,  $\omega_{\text{res}}$  and  $\omega_{\text{cut}}$  depend on  $R$  schematically as shown in Fig. 2. For a given frequency,  $\omega$ , the position of the resonance and cut-off,  $R_{\text{res}}$  and  $R_{\text{cut}}$ , are uniquely determined by  $f$ .

The index of refraction profile is like that shown in Fig. 3. The fast wave launched inward from the outer boundary is reflected at  $R = R_{\text{cut}}$  and goes back to the outer boundary. A receiver B loop

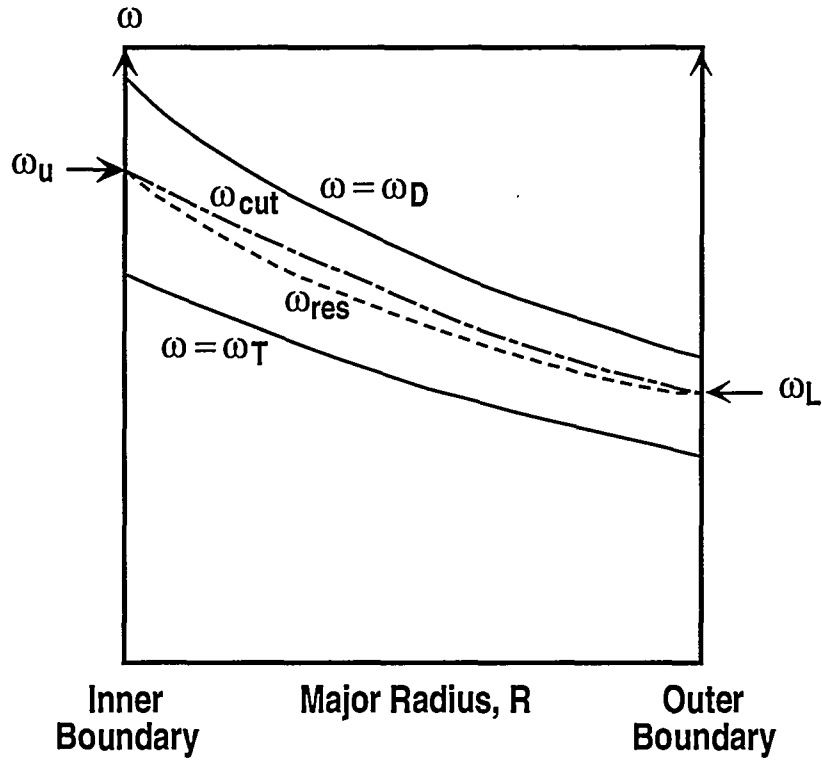


Fig. 2. The ion cyclotron frequencies of deuterium and tritium,  $\omega_D$  and  $\omega_T$ , and resonance and cut-off frequencies as a function of tokamak major radius.

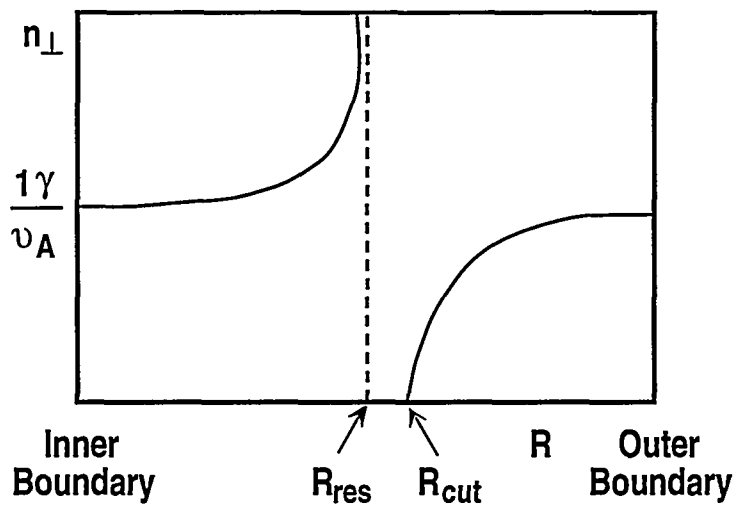


Fig. 3. The Alfvén index of refraction,  $n_{\perp} \equiv v_A k_{\perp} / \omega$ , as a function of major radius.

outside the outer boundary picks up the reflected wave. The measurement of this round trip phase shift (reflectometry) finds  $R_{\text{cut}}$  with the aid of the mass density information obtained by the previously described diagnostic.

By sweeping the frequency, the profile of  $f$  may be obtained. Figure 2 also indicates that the wave propagates across the plasma at the frequencies  $\omega < \omega_L$  and  $\omega > \omega_U$ . The lower end of the absorption band  $\omega_L$  gives  $f$  at the outer boundary and the upper end frequency  $\omega_U$  determines  $f$  at the inner plasma boundary.

We have made the computer code to invert the phase vs. frequency data to the species mix as a function of radius. We have adopted a simple numerical method given by Doyle, et al.<sup>3</sup> The computation time is short enough to process the data between the usual tokamak shots.

The proof-of-principle experiments of the reflectometer will be carried out on the DIII-D tokamak in the near future.



## REFERENCES

- <sup>1</sup>H. Ikezi, R. I. Pinsky, S. C. Chiu, and J. S. deGrassie, "Fast Wave Propagation Studies in DIII-D Tokamak," accepted for publication in *Phys. Plasma*.
- <sup>2</sup>T. H. Stix, *Waves in Plasmas*, American Institute of Physics, New York (1992) p. 7.
- <sup>3</sup>E. J. Doyle, T. H. Lehecka, N. C. Luhmann, Jr., W. A. Peebles, and the DIII-D Group, *Rev. Sci. Instrum.* **61**, 2896 (1990).

